


The Bobath is Esthetica's medical treatment table, meticulously engineered for safety, stability, and patient comfort in neurological rehabilitation and therapy environments. Its sturdy 2-section design features an elongated backrest outfitted with two premium gas springs, ensuring effortless and comfortable patient positioning throughout every therapy session. Available in two variants — the Pro model, equipped with swivel wheels and a walk-free switch for seamless navigation in busy clinical environments, and the Basic model, featuring a foot switch for straightforward, stationary operation — the Bobath delivers reliability, safety, and clinical versatility across both configurations.



KEY FEATURES

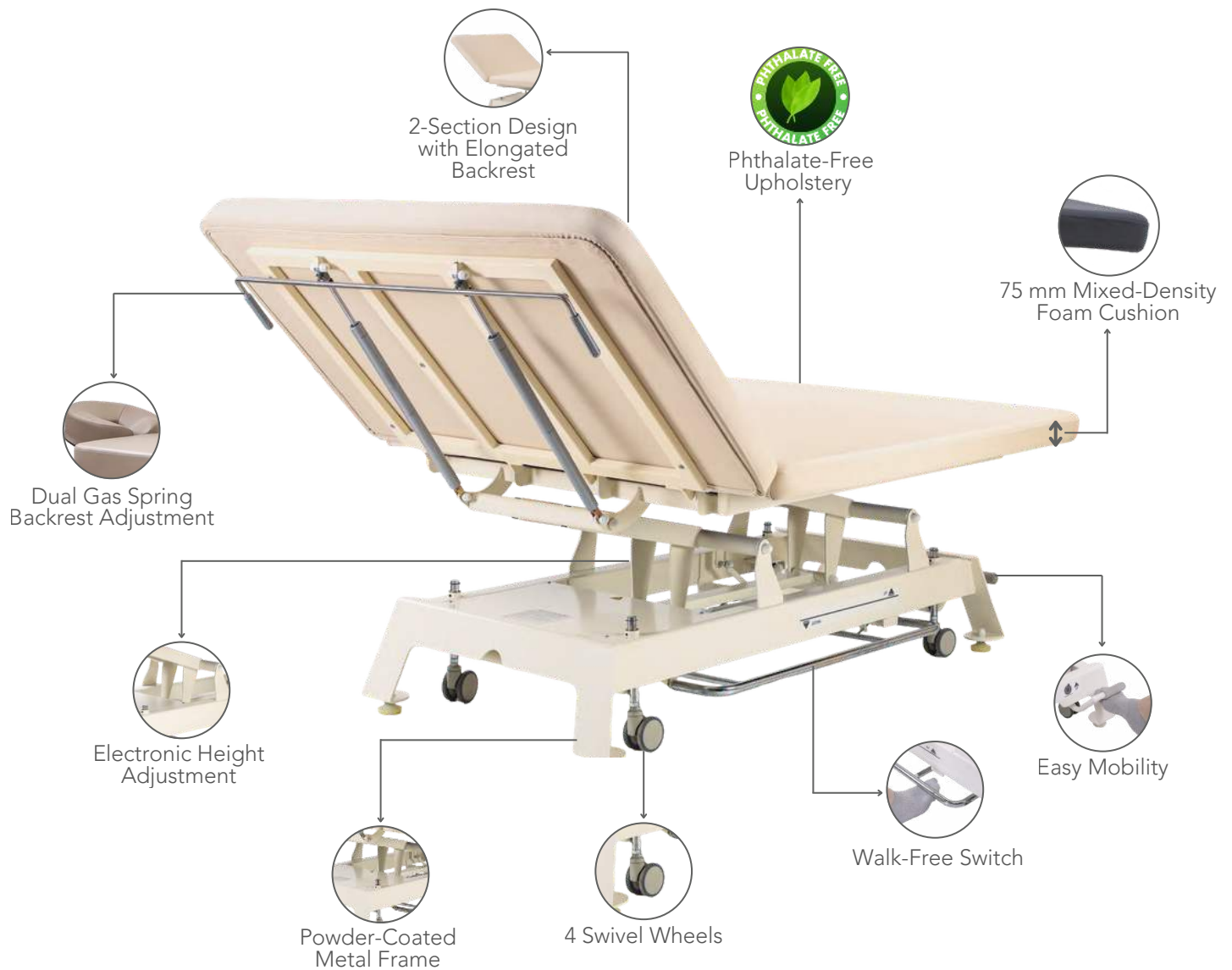

2-Section Design with Dual Gas Spring Backrest

The Bobath's 2-section design features an elongated backrest fitted with two premium gas springs — delivering effortless, comfortable patient positioning that supports optimal therapy outcomes across neurological rehabilitation sessions.

Phthalate-Free Upholstery

The Bobath is upholstered in phthalate-free fabric designed for clinical environments — safe, durable, and built for the demands of daily professional use, with antibacterial, fire retardant, and stain resistant properties as standard.



Cream



Beige



Brown



Mouse Grey

Electronic Height Adjustment

Smooth electronic height adjustment allows clinicians to set the ideal working height quickly and consistently across every session — reducing physical strain and maintaining ergonomic comfort throughout the treatment day.

Walk-Free Switch & Swivel Wheels

The walk-free switch allows the Bobath to be operated from both sides of the table without the clinician needing to reposition, while swivel wheels allow seamless navigation through busy clinical environments — streamlining workflows and maximising treatment efficiency.

FEATURES

Design	Two piece articulating top
Backrest	Manually Adjustable (Gas Spring)
Height	Electronically Adjustable
Walk Free Switch	Enables you to operate the table from both sides
Castor Wheels	For easy mobility in the treatment room

TECHNICAL SPECIFICATIONS

Length	2040 mm
Width	1220 mm
Minimum Height	530 mm
Maximum Height	890 mm
Backrest Inclination	-25° (when lowered) & 96° (when raised)

ELECTRICAL PARAMETERS

Power in (Mains Power Supply)	100-240 V, 50-60 HZ, 2A
IPX Rating	IPX4
Operating Temperature	5-40°C
Shock Protection:	Class I

BASE & FRAME

Frame Structure	Mild steel construction powder coated frame & ABS molded parts
Base	Powder coated metal frame
Safe Working Load	180 kgs

SAFETY STANDARDS

General safety for medical electrical equipment	IEC 60601-1
EMI/EMC	IEC 60601-2

UPHOLSTERY


Type	Phthalate-Free Upholstery
Thickness	1.01 mm
Tensile Strength	ISO 1421, 804 N/ 5 cm
Colour Options	Cream, Beige, Brown, Mouse Grey
Weight	732 g / sqmt
Abrasion	ISO 5470-2 = 2,00,000 cycles / min

WARRANTY


3 years on motors





TECHNICAL DRAWINGS

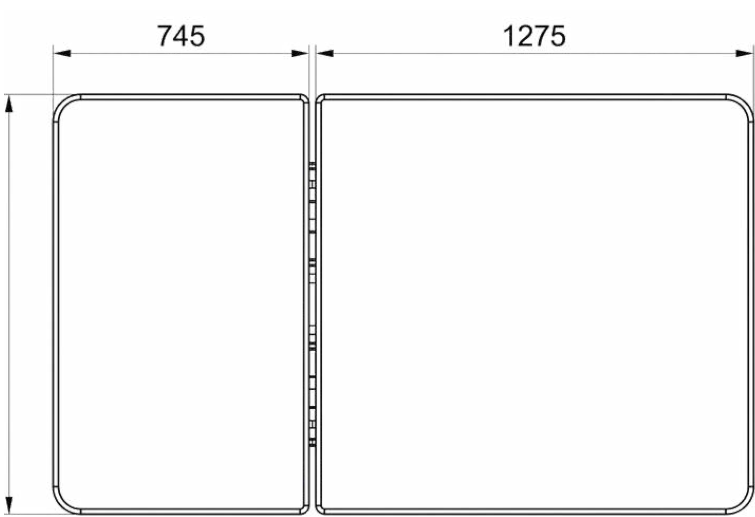


TOP VIEW

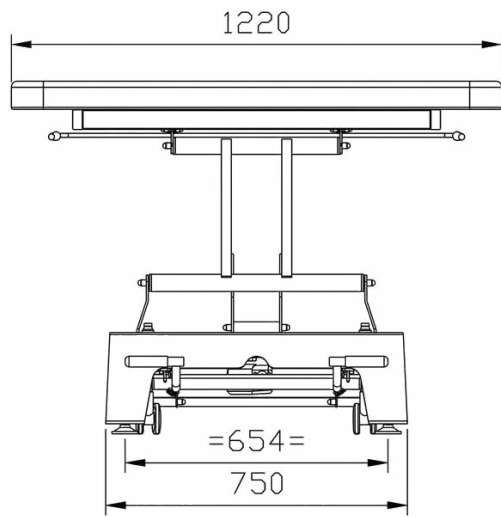


FRONT VIEW

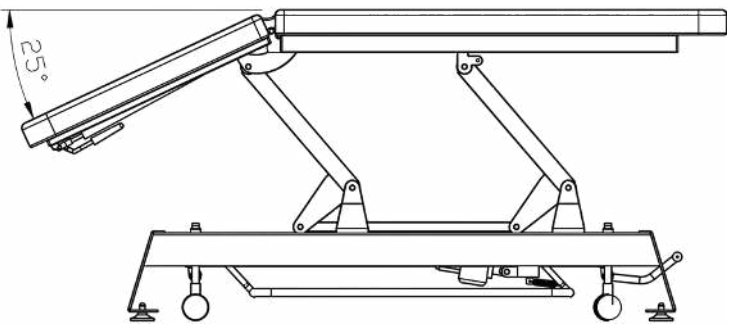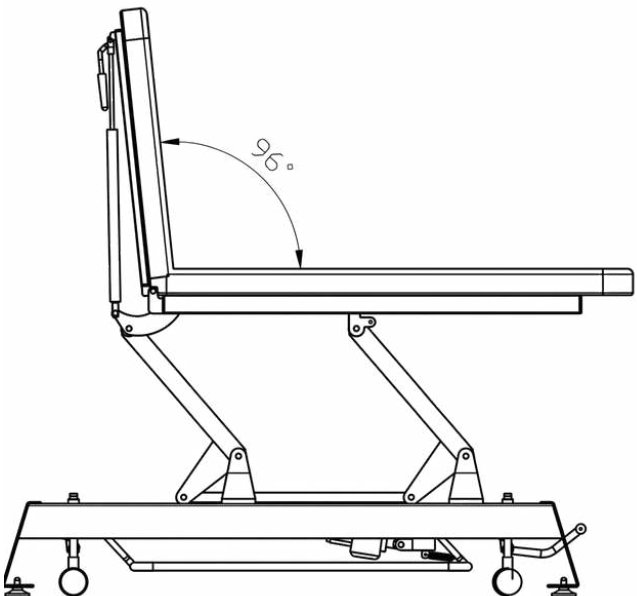

BACKREST INCLINATION

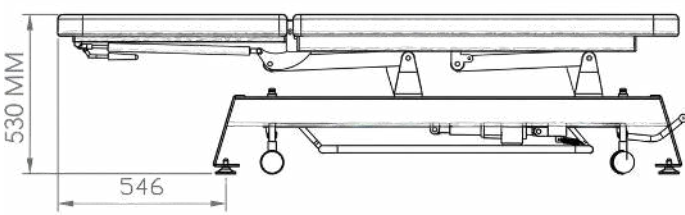

MINIMUM HEIGHT

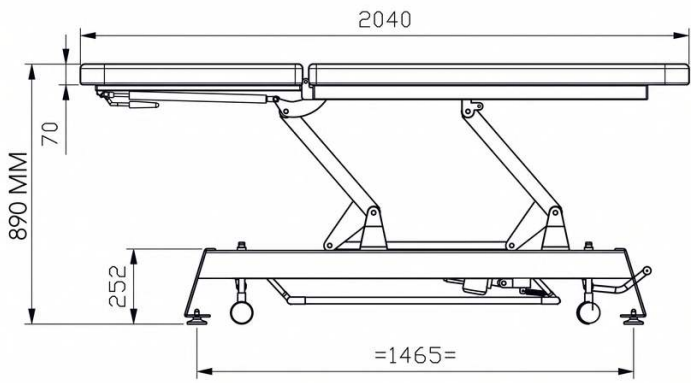

MAXIMUM HEIGHT

ALL DIMENSIONS ARE IN MM