

ESTHETICA BASIC PHLEBOTOMY CHAIR

The Esthetica Basic is a reliable, ergonomically designed stationary phlebotomy chair for hospitals, clinics, diagnostic laboratories, and blood banks. A height-adjustable blood drawing arm provides proper positioning for safe and accurate venous access, while a padded flip armrest with height and depth adjustment and an adjustable elbow rest ensure patient comfort and stability throughout every procedure. The flip armrest swings away smoothly for easy patient exit. Built on a heavy-duty mild steel sheet metal and tube frame with easy-clean phthalate-free upholstery, the Basic delivers the durability and hygiene that high-volume clinical use demands.



KEY FEATURES



Height Adjustable Blood Drawing Arm

On the draw side, the dedicated blood drawing arm adjusts in height to bring the patient's arm to the ideal level for venous access, giving the phlebotomist a clean, stable working surface for blood draws and other phlebotomy procedures.

Phthalate-Free Upholstery

The Esthetica Basic is upholstered in phthalate-free upholstery, designed for professional clinical environments. It resists stains, oils, and disinfectants, withstands frequent cleaning, and carries antibacterial and fire retardant properties as standard.



Padded Flip Armrest with Height & Depth Adjustment

On the support side, a padded flip armrest holds the patient's resting arm in a comfortable, secure position. Adjustable in both height and depth to suit each patient, it swings clear of the seat for easy, unassisted entry and exit.

Heavy-Duty Mild Steel Frame

Constructed from mild steel sheet metal and tube, fully welded and powder-coated, the frame provides a stable, rigid platform built to withstand years of continuous clinical use without flex or wear.

FEATURES

Armrests	Blood Drawing Armrest & Padded Flip Armrest with Height & Depth Adjustment
Backrest	Padded Backrest
Footrest	Fixed
Height	Fixed-Height Stationary Base

TECHNICAL SPECIFICATIONS

Backrest Dimention	490 mm (L) × 360 mm (W)
Seat Dimension	480 mm (L) × 530 mm (W)
Padded Flip Armrest	620 mm (L) X 150 mm (W)
Blood Drawing Armrest	435 mm (L) X 140 mm (W)
Seat Height	530 mm
Overall Height	910 mm

BASE & FRAME

Base	Heavy-Duty Mild Steel Sheet Metal & Tube Frame
Safe Working Load	200 kgs
Patient Load Bearing capacity	150 kgs


UPHOLSTERY

Type	Phthalate-Free Upholstery
Thickness	1.01 mm
Tensile Strength	ISO 1421, 804 N/ 5 cm
Colour Options	Cream, Beige, Brown, Mouse Grey
Weight	732 g / sqmt
Abrasion	ISO 5470-2 = 2,00,000 cycles / min


WARRANTY

1 years on structure


TECHNICAL DRAWINGS

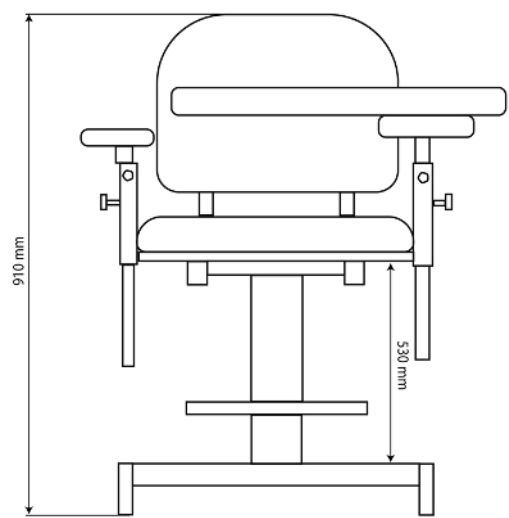


FRONT VIEW

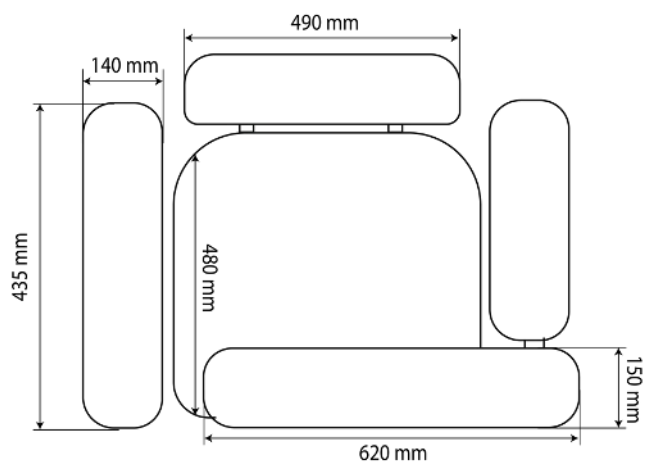

TOP VIEW

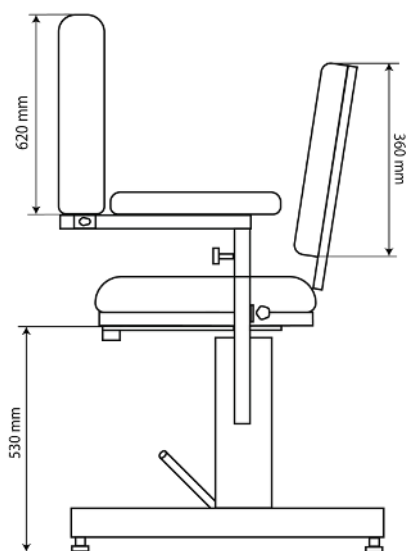

SIDE VIEW

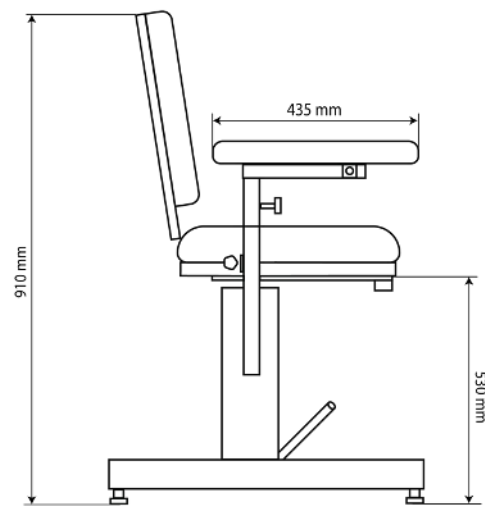

HEIGHT

ALL DIMENSIONS ARE IN MM